


adnar
by D:O IMPLANT

Dental Treatment Unit

ADNAR Prime

The chair comes in a hanging hose dentist element configuration with push button control. It has double armrest, advanced foot control, 8 bulb LED operating light and an inbuilt water heater. The chair comes with our standard PU upholstery.

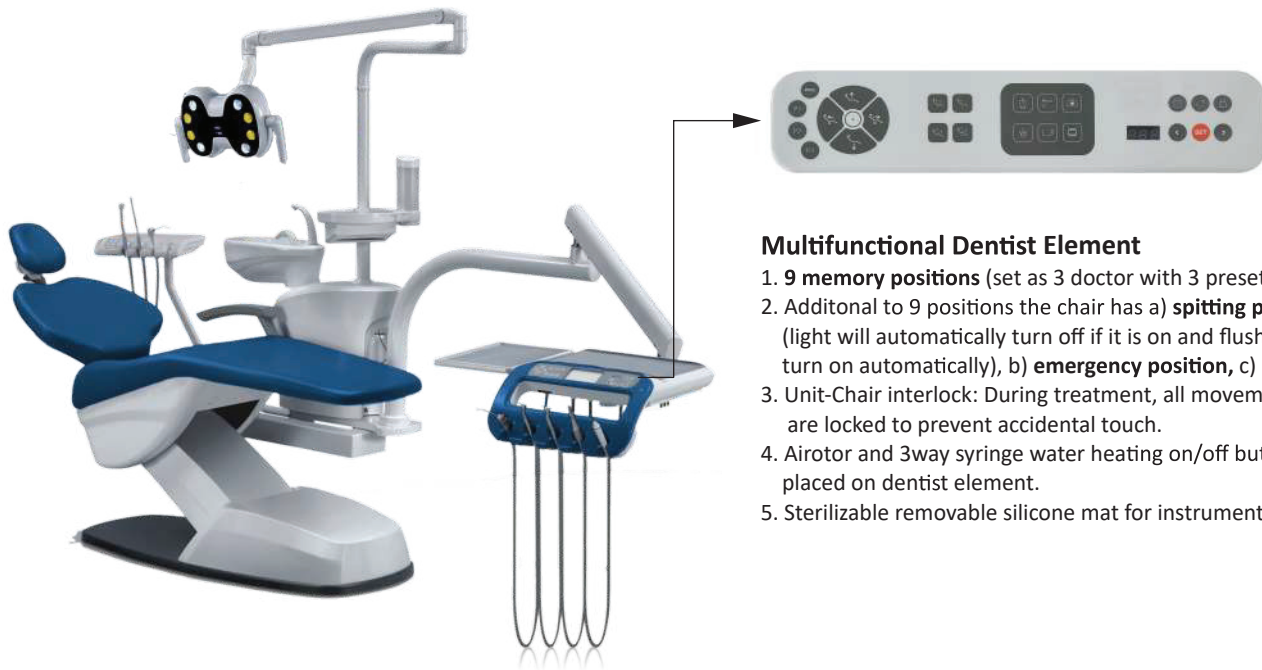

Multifunctional Dentist Element

1. **9 memory positions** (set as 3 doctor with 3 presets)
2. Additional to 9 positions the chair has a) **spitting position** (light will automatically turn off if it is on and flushing water will turn on automatically), b) **emergency position**, c) **return to last position**
3. Unit-Chair interlock: During treatment, all movements of the chair are locked to prevent accidental touch.
4. Airtor and 3way syringe water heating on/off button is conveniently placed on dentist element.
5. Sterilizable removable silicone mat for instrument.



8 Bulb LED
Operating Light



Detachable
Rotary Spittoon



Low Voltage Heater



Assistant Control



3 in 1 Holder



Advance Foot Pedal

ADNAR Progress

The chair comes in a hanging hose dentist element configuration with touch button control. It has a 12 bulb LED operating light, stainless steel instrument tray, automatic disinfection system, wireless phone charger, hidden face mirror, and an inbuilt water heater. The chair comes with floor type water unit and has a soft seamless leather upholstery

ADNAR Pinnacle

The chair comes in a hanging hose dentist element configuration with 12.3" touchscreen for controlling the chair positions, micromotor and scaler setting and record reports. It has a 12 bulb LED operating light, stainless steel instrument tray, automatic disinfection system, wireless phone charger, hidden face mirror, water sterilizer for tap water and an inbuilt water heater. The chair comes with floor type water unit and has a soft seamless leather upholstery



Wireless charging



12 Bulb LED Operating Light



Faro EDI LED Light (Optional)



HH Dentist Control



HH Dentist Control




Assistant Control




Floor Mount Water Unit



Water Sterilizer (Optional)

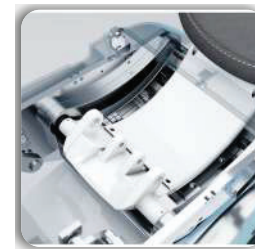

Soft Movement Synchronized

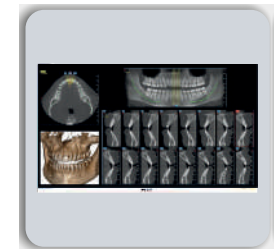

Monitor with Arm (Optional)

ONE-CLICK INTELLIGENT DISINFECTION

Available in Adnar Progress and Pinnacle

Digital automatic waterway disinfection system

- Fully inclusive: covers all water pipelines for diagnosis and treatment to ensure comprehensive disinfection.
- Fully automatic: one-click fully automatic disinfection, no need for on-duty, efficient and worry-free.
- Fully intelligent: visualize the disinfection process and status, accurately control, optimize the process.
- Full record: disinfection records can be viewed, saved up to 5 years, exported, and deleted.
- Fully up to standard: meets the national requirements for dental unit pipelines and treatment water.



COMFORTABLE EXPERIENCE FOR DENTISTS AND PATIENTS





Dental Treatment Unit



D:O IMPLANT
DIO DIGITAL IMPLANT INDIA PVT. LTD.

Gurgaon Office: 7061, 7062, 708,
709, 7101, 7102, 7111 & 7112,
7th Floor, Block B, Vipul Square,
Sushant Lok, Phase I,
Gurgaon, Haryana - 122009
Contact: 0124-4751803,
E-mail: info@dioimplant.co.in

Mumbai Office: Unit No. 605,
Business Suites 9, Plot No. 93
S.V. Road, Santacruz,
West Mumbai - 400054
Contact: 022-46021803

Bangalore Office: No.18/6, First Floor,
Cunningham Road,
Bangalore - 560052

Chennai Office: No. 45,
3rd Floor Ashok Towers,
KB Dasan Road, Nandhiniee Sweets,
Alwarpet, Chennai - 600018

Website: www.dioimplant.co.in